
SEAGO

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
Thanks for choosing SEAGO.
Before use, please read the instructions carefully
and keep it at a safe place for future reference.

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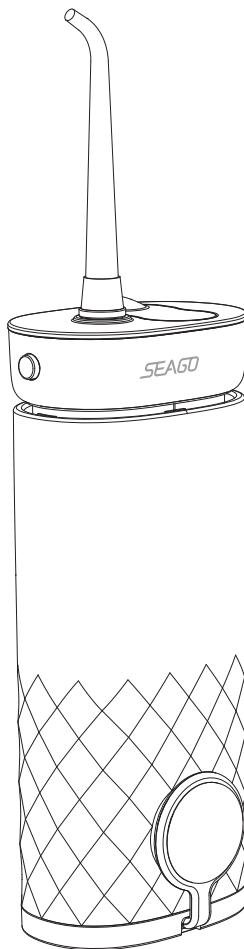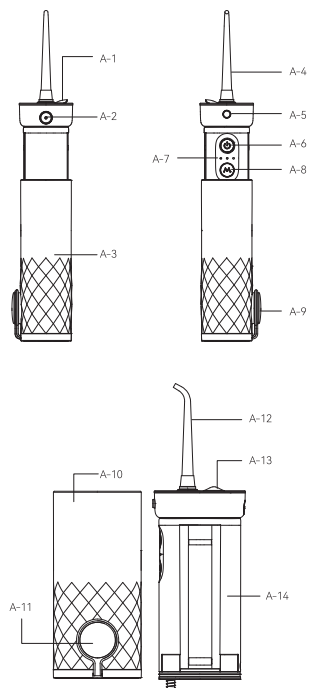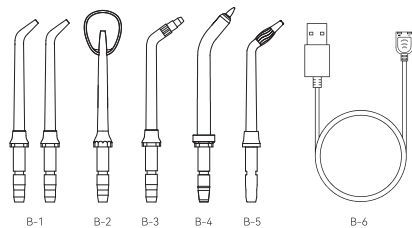
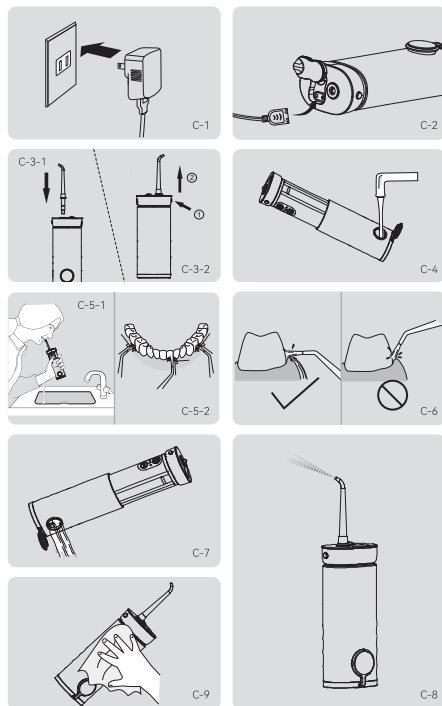
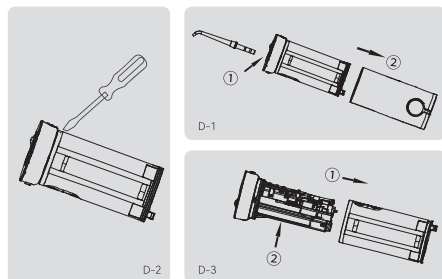
Safety Instructions

Warning

1. Do not use wet hands to connect the charging cable interface or unplug the power adapter. And do not use wet hands to connect or unplug the power adapter from the household power outlet.
2. Please make sure the charging port is dry when charging! Please use the original charging cable!
3. The product is not designed for someone with physical, sensory or mental disabilities or lack of experience and knowledge such as children. They can use it with the help of guardian. And guardian is responsible for the safety of users and take care of children and make sure they don't play with the product.
4. This product contains a battery, do not put it into the fire to heat, or charge, use or place in a high temperature place, and avoid direct sunlight.
5. Do not immerse the product in water for a long time to avoid product function failure.
6. Do not expect damage, damage threat to modify the charging cable. Don't place heavy objects on the cable or squeeze the cable during charging.
7. Do not knock or drop the product, if this product is damaged due to collision and fall, please stop using it;
8. If there are any abnormality and failure, please stop using or charging the product.
9. When cleaning the power adapter and charging cable, please be sure to remove the power adapter from the household power outlet.
10. Please clean the charging line interface and power adapter plug regularly every 6 months to prevent dust accumulation.

Attention

1. Do not use the product without installing a nozzle.
2. This product should not be used by people with no operation-ability. This product should not be used by people who are not sensitive to oral sensation.
3. Please don't put the product where child or baby can reach. And, please prevent them touching or operating the product.
4. Patients with periodontal disease or those who have undergone oral surgery in the past 2 months, please consult a doctor before use.
5. Do not use the product for any purpose other than oral cleaning.
6. Do not fill the water tank with water over 40°C, otherwise it may cause burns.
7. Do not share the nozzle with families or others.

8. Only the nozzles and accessories provided, sold and recommended by Seago are allowed.
9. Please do not wrap the charging cable on the body when storing.
10. Please do not pull the charging cable and hold the power adapter and pull it out when disconnecting the charging plug.
11. Please clean and maintain the equipment after use. (reference routine maintenance)
12. Please do not repair and maintain it when using the product.
13. If the power plug or power cord is damaged, please replace it.
14. Please don't refit, disassemble or repair the product by yourself.
15. Do not fill the product with iodine, saline, or concentrated essential oils that are insoluble in water. Using these liquids will reduce product performance and shorten service life; please do not fill pure mouthwash liquid into this product, if you need to use mouthwash, please dilute it in proportion and then fill it.

Component Name

WATER FLOSSER

- A-1 Magnetic Charging Port
- Water Injection Port
- Water Tank Cover
- A-2 Pause Button
- A-3 Retractable Water Tank

A-4 Nozzle

- A-5 Nozzle Release Button
- A-6 Switch / Charge Indicator
- A-7 Mode Indicator
- A-8 Mode Button
- A-9 Water Injection Port

DETACHABLE PART

- A-10 The Water Tank
- A-11 The Tank Cover

A-12 Nozzle

- A-13 The Tank Cover
- A-14 Main Body

NOZZLE

- B-1 Standard Nozzle:** for daily clean
- B-2 Tongue Scraper:** for tongue cleanliness and freshness
- B-3 Orthodontic Nozzle:** for individuals with false teeth, dental braces, and dental implants
- B-4 Periodontal Sprinkler:** for periodontal pocket cleaning
- B-5 Tornado Nozzle:** Rotary friction cleans residues and plaque
- B-6 Magnetic Charging Cable**

Product Parameters

Material	ABS(body), PMMA(nozzle)
Body Size	About 151×47×63mm
Body Weight	About 230g
Rechargeable Battery	Lithium-ion
Charge Mode	Magnetic charging
	5-40°C (under standard atmospheric pressure)
Applicable Altitude	Below 5000m
Charging Time	About 4 hours
Nozzle	6pcs
Mode	3 Mode : DIY / Soft / Clean
Hours of Continuous Use	75 times (Clean mode, one box of water at a time)
Power Consumption	2W
Input Voltage	5V
	Under 90%
Waterproof Level	IPX7

	Standby
	Insulation type: Class BF equipment
	Remind readers
	Upper and lower temperature limits
	Humidity upper and lower limits
	EU environmental protection directive
	Manufacturer name and manufacturer address
	Production batch
IPX7	Waterproof level
	DC symbol

Use Guide

• Before use: please follow the steps below to charge the product:

1. Product application part is BF type.

1. Check and make sure the product and hands are dry, then connect the magnetic charging cable to the power supply adapter and plug the power plug into the power supply socket. (C-1)
2. To avoid toppling during charging, you need to push the water box back to the oral irrigator and horizontally place it for charging. The charger plug should be put into the charge port of the oral irrigator properly (C-2).
3. An adapter of 5V, 1A (actual maximum output power less than 15W) that conforms to GB 4705.1 – 2005 shall be used to recharge the device.

1. TIPS

When the red indicator light flashes or the product cannot be turned on normally after use, it means that the battery power is low and the product needs to be charged.

- It takes about 4 hours to fully charge.
- When charging, the red charging indicator light is on.
- The green charging indicator light is on, indicating that the charging is complete.
- When the battery is fully charged, the charging will be turned off for battery overcharge protection.
- During charging, the red indicator flashes quickly, indicating that the battery is overheated. Please stop charging.

When using, please follow the steps:

1. Insert the nozzle into the hole at the top of the handle (C-3-1). If the nozzle needs to be replaced, shut down the water flosser first and then press the release button of the nozzle lightly to easily replace the nozzle (C-3-2).
2. Open the water tank cover and fill the water tank (C-4). If it needs to increase the water storage capacity, it can open the water tank cover and then pull the water tank. Note that after the water tank is pulled to the stop line, it can't continue pulling.
3. This product has three modes: DIY / Soft / Clean.
First Mode (DIY mode): Continuously sprays water, the first indicator light on the left lights up, and the timing is 2 minutes. After the first mode is turned on, press and hold the gear button for a long time, and the water spray intensity will change circularly from weak to strong, with a cycle period of 5 seconds. After releasing the button, the setting intensity will be memorized.
Second mode (Soft mode): Continuously sprays water, two indicator lights on the left light up, and the timing is 2 minutes.
Third mode (Clean mode): Continuously sprays water, three indicator lights light up, and the timing is 2 minutes.

Button instructions:

1. Switch button: Press the switch button to turn on the product, and press the switch button again to turn off the product.
2. Mode button: short press the mode button to shift modes in the shut down and start up state (shift mode in the shut down state, and the product remains non-operating). The shift order is first-second-third-first.....
3. Spot-spray function: When the product is in working condition, press the pause button to stop spraying water, and release it to resume spraying water.

1. TIPS

If you use our product for the first time please try DIY mode for a week until adaptation. Then try other modes.

4. Before use, hold the water flosser vertically, with the nozzle sticking into the mouth and aiming at the teeth or gums. At the same time, open the mouth slightly so that water can flow out smoothly (C-5-1&2).
5. Press the switch button to use.

Please note the following points when using:

- A) align the gum line with the nozzle and be perpendicular to the gum;
- B) move slowly along the teeth and clean them in sequence;
- C) clean teeth, trays, crowns and dentures thoroughly.

① TIPS

Do not spray water directly down the periodontal pouch (C-6).

After use, please press the switch button to turn off the machine.

This product has 2 minutes timing function, use 2 minutes automatic shut down. If you still need to use it, please press the switch button.

The product has the mode memory function.

After use, please follow the following steps:

1. Open the tank cap and pour the water out (C-7).
2. Press the switch button, the water flosser works, and automatically drain the water inside (C-8).
3. Press the switch button to turn off the water flosser and then close the tank cover.
4. Dry the water flosser with a cloth (C-9).

① TIPS

To avoid the spread of bacteria, make sure there is no residual water in the water flosser.

If you plan not to use the product for a long time, please stock the product after dry.

Daily Maintenance

1. Please clean the product before using. Or cleaning the product every 6 months if you don't use it for a long time.

2. Clean the product with clean water or neutral detergent. Do not use corrosive or abrasive detergent. Do not use hot water over 40°C in order to avoid damaging the product.

Cleaning the body:

- ① Please clean the body with dry or wet cloth.
- ② Do not wash the body into the water.

Cleaning of water tank:

- ① Wipe off with a dry or wet cloth, or remove and rinse with clean water.
- ② If you plan not to use it for a long time, please clear the water.

Cleaning of nozzle:

- ① You can use a dry cloth or a wet cloth to wipe, or take it off and rinse with water, then dry with a cloth.
- ② Do not bend, pull and twist nozzle pipes.
- ③ Dental experts recommend. Replace the nozzle every 6 months.

Trouble Shooting

The Problem	Probable Cause	The Solution
Can't work	Newly purchased or not used for more than three months	Please charge for more than 4 hours
The water pressure is low at work	Nozzle blocked by something Nozzle deformation	Replace the nozzle
	The screen is blocked	Clean the screen
It can only be used for a short time after full charge	The battery is at the end of its useful life	Please contact customer service
	Insufficient charging time	Please charge for more than 4 hours before use
The charging indicator does not work when charging	The power connection is not properly connected	Please connect the power correctly
	There is a fault in the charging circuit of the product	Please contact customer service
No water overflowing	New purchases or unused for a while	Pull out the nozzle to flush the nozzle hole on the main body

Scrap Processing

Remove the rechargeable battery

Do not disassemble the product by yourself unless you want to dispose of it premature removal will damage the product. Operation is subject to basic safety precautions. Before removing the battery, you must unplug the power plug and charging plug.

1. Before removing the battery, make sure the battery is exhausted. If you need to run out of power, you can press the switch button to turn on the water flosser and let it work until it stops automatically.
2. Press the nozzle release button.(D-1-①), remove the nozzle and forcibly remove the water tank in the direction of the arrow (D-1-②).
3. Use a flat screw driver to pry it open from (D-2).
4. Take off the shell as shown (D-3-①), and remove the battery (D-3-②).

⚠ Attention

- ① Properly recycle or dispose of the battery. Ask your local waste management department whether to recycle other parts of the product.
- ② Be careful not to short the positive and negative terminals of the battery. Don't charge the product after removing the battery.

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Dear customers, thank you very much for your concern and support for our electric toothbrush, water flosser and other products! Thank you for choosing and purchasing our products! The company promises that this product has passed our company's strict quality inspection. If there is any quality problem from the date of purchase, it will be returned within 7 days, replaced within 15 days and guaranteed within 2 years. If the operation is not in accordance with the instructions or without the company's authorization units and individuals to open or repair, is not included in this scope.

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Fill out information below by customer

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Fill out information below by retailer

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