XF7小矮
SEAGOO 碧碧

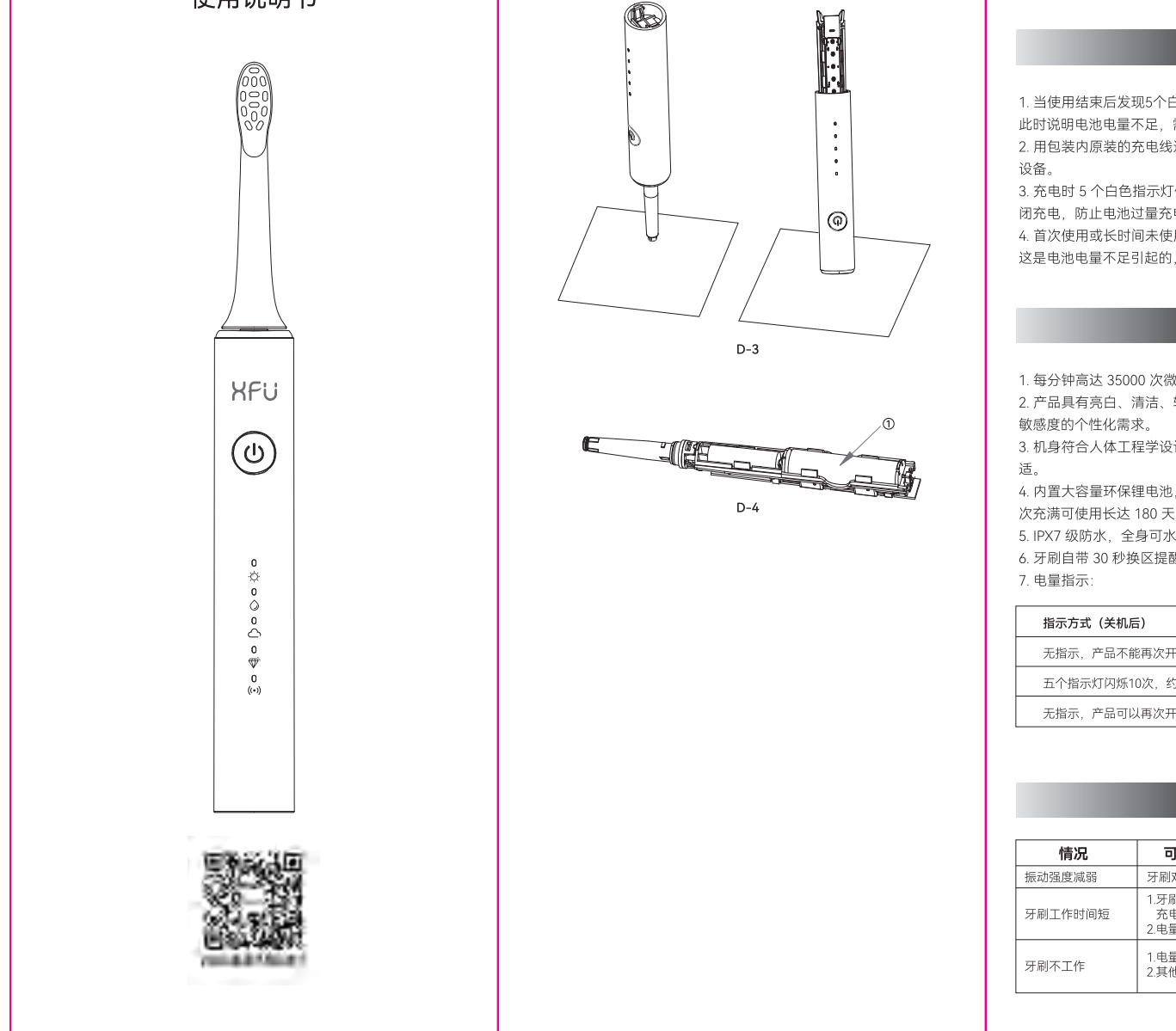
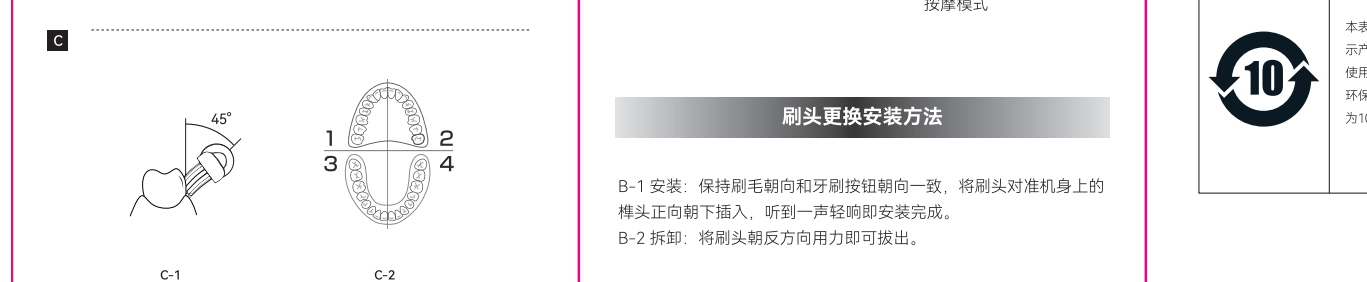
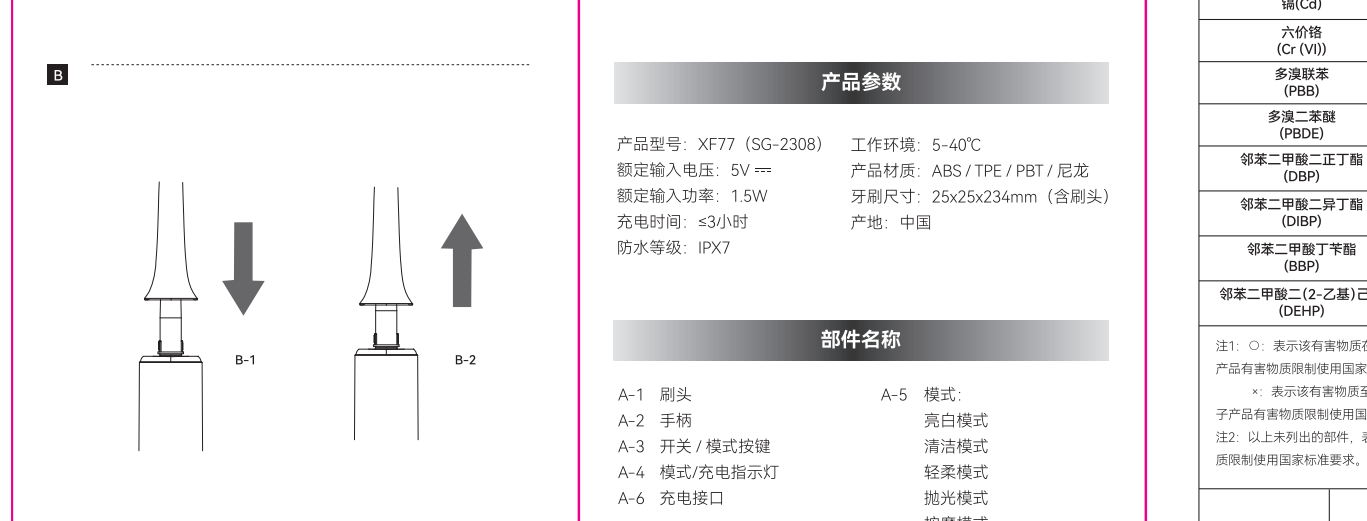
声波电动牙刷
Electric Toothbrush
XF77 (SG-2308)

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Instruction Manual
使用说明书

[illegible][illegible]

1. 启动前，请仔细阅读本说明书，并严格按照本说明书的要求进行操作。 2. 使用过程中，请勿将本产品置于潮湿、高温、阳光直射等环境中。 3. 请勿将本产品用于非设计用途，以免造成人身伤害或财产损失。 4. 请勿将本产品用于易燃易爆场所，以免引发火灾或爆炸。 5. 请勿将本产品用于医疗、工业、科研等特殊用途，除非另有说明。 6. 请勿将本产品用于非法用途，以免承担法律责任。 7. 请勿将本产品用于高风险作业，以免发生安全事故。 8. 请勿将本产品用于超出其额定功率的场合，以免损坏产品。 9. 请勿将本产品用于超出其额定电压的场合，以免损坏产品。 10. 请勿将本产品用于超出其额定温度的场合，以免损坏产品。	1. Turn off the appliance when leaving it unattended. 2. Do not use the product near flammable or explosive materials or in the same environment. 3. If you have undergone oral or dental surgery, use the product 2 weeks after the procedure to avoid bleeding. 4. The product is intended for general cleaning, gums, and the tongue only. Do not use it for any other purposes. 5. Do not use the product with deformed bristles. Replace the brush head immediately if signs of wear appear. It is recommended to replace the brush head every three months. 6. Use only genuine brush heads produced by our company. If you are using incompatible accessories, such as sodium disodium bicarbonate, peroxide, or other carbonate ingredients (common whitening toothpastes), thoroughly clean the brush head after use to prevent plastic crystallization.	1. Turn off the appliance when leaving it unattended. 2. Do not use the product near flammable or explosive materials or in the same environment. 3. If you have undergone oral or dental surgery, use the product 2 weeks after the procedure to avoid bleeding. 4. The product is intended for general cleaning, gums, and the tongue only. Do not use it for any other purposes. 5. Do not use the product with deformed bristles. Replace the brush head immediately if signs of wear appear. It is recommended to replace the brush head every three months. 6. Use only genuine brush heads produced by our company. If you are using incompatible accessories, such as sodium disodium bicarbonate, peroxide, or other carbonate ingredients (common whitening toothpastes), thoroughly clean the brush head after use to prevent plastic crystallization.
电容量约 100% 电容量约 70%~90% 电容量约 70%~90%	Model: XF7T (SG-238) Rated Input Voltage: 5V ~ 30V Rated Power: 1.5 W Operating Time: 4~3 hours Charging Temperature: 5~40°C (41~104°F) Product Dimensions: 25 × 28 × 28mm (including brush head)	Product Specifications Waterproof Rating: IPX7 Material: ABS+TPC (PET) / TPU Country of Origin: China
产品参数	快速指南 1. 准备工作 2. 充电 3. 使用 4. 清洁 5. 存储	Parts Name A-1: Brush Head A-2: Handle A-3: Mode / Power Button A-4: Mode / Charging Indicator Light A-5: Charging Port A-6: Modes A-7: Whiteness Mode A-8: Clean Mode A-9: Sensitive Mode A-10: Push Mode A-11: Massage Mode

Indicator	Behavior	Battery level
	No indication product cannot be turned on again	Less than 10%
	Five indicator lights flash 10 times (usually, 10 seconds)	Approx. 10% - 30%
	No indication, product can be turned on normally again	Approx. 30% - 100%

Trouble Shooting		
Issue	Possible cause	Solution
Reduced vibration intensity	Excessive pressure applied to the tooth	Gently push the bridge against your tooth without pressing hard
Shortened operating time	1. Impinger pollution leading to poor battery performance 2. Insufficient battery level	1. Recharge the toothbrush completely on a regular basis 2. Extend the charging time
Device does not turn on	1. Battery is completely drained 2. Other internal issues	1. Charge the device 2. Contact SEAOI customer support

Battery Disposal

This product uses a lithium battery. When the product reaches the end of its life, you must remove the battery and dispose of it properly according to local waste management regulations. Please contact your local waste disposal authority for specific recycling instructions.

Tools required for battery removal: Saw, Wire Cutters.

Warning: Only remove the battery when discarding the product. Before disassembly, ensure the device is disconnected from the power supply and turned off. Completely drain the battery before removal. When using tools, strictly follow basic safety rules and take appropriate protective measures to avoid injury.

1. Remove the brush head. Use a saw to cut through the handle at the marked position (20mm from the top) (D-1) (D-2) (Battery).
2. Image shown: Sawing effect / cross-section (D-1) (D-2).
3. Place the connecting post of the cross-section against a flat surface (such as a table), then hold the handle and press down firmly to separate the brush head from the handle.
4. Once the bracket is removed, cut the battery wires and take out the battery (D-4).

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Bruah Head Installation and Removal

1. Installation: Align the bridge direction with the tooth direction on the handle. If the bruah point goes onto the wronging spot of the first main body and press it straight down until you hear a click, indicating it is securely installed.

2. Removal: Pull the bruah upward with force to detach it from the handle.

Operating Instructions

1. Install the bruah firmly onto the handle. Apply a suitable amount of toothpaste onto the bristles.

2. Hold the handle at the base of your teeth and gums, holding the bruah at a 45°angle to the tooth surface (see Fig. 1). Press the power button to start. Slightly move the bruah horizontally across the tooth surface.

3. Divide your mouth into four sections: lower left front, upper right, upper left, lower right. Spend 10 seconds brushing each side of the inner and chewing surfaces of each section (see Fig. 2).

4. The toothbrush will turn around automatically to work from the top to the next quadrant. The appliance will automatically stop after 2 minutes of continuous use.

5. Within 3 seconds of turning on the device, pressing the button will switch the mode. The switching sequence is: Whiskering → Sensation → Massage → Rotate. When the mode is changed, if the mode is pressed within three seconds, the button stage will skip the off delay.

6. The mode used in the current session will be memorized and the same mode will be used the next time you turn on the device.

7. Rinse the bruah thoroughly with clean water after each use. Do not soak, but wash the handle both the diverging gear "light-colored" waterside and "dark-colored" sides for a deep clean.

- 1. When the 5 white indicator lights on the product cannot be turned on normally after use, the battery is low and needs to be charged.
- 2. Connect the toothbrush to a USB charger or a device with a USB port using the original charging cable provided in the package.
- 3. During charging, the 5 white indicator lights will illuminate one after one in sequence. When fully charged, all lights will turn off and charging will stop automatically.
- 4. If the indicator light does not illuminate during the first 5 minutes of charging (especially during the first 2 periods of 1.5 periods each), this is a normal self-protection mechanism due to relatively low battery level. The lights will turn on automatically after 5 minutes.

Features

- 1. Up to 1300 maximum vibrations per minute for a deep and thorough cleaning.
- 2. 27 timing modes (Whitening, Clean, Sensitive, Polish, and Massage) to meet personalized gum sensitivity needs.
- 3. Ergonomic design that is perfect for the hand and single-hand operation for a comfortable grip.
- 4. 1500mAh capacity, one-year (12 months) battery life and single-charge lasting up to 180 days (based on 27 minutes of use, twice a day as 15-second mode).
- 5. IPX7 waterproof rating allows for safe use in the shower and easy cleaning under running water.
- 6. Smart timer with a 30-second quadrant start and a 2-minute auto-shutdown plus a convenient memory function.
- 7. Battery Level Indicator (after charge stop)

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Dear customers, thank you very much for your concern and support for our electric toothbrush, water flosser and other products. Thank you for choosing the products of the company. The company promises that the water flosser has passed our company's strict quality management system and will be produced from the date of purchase, it will be returned within 7 days, replaced within 15 days and guaranteed within 2 years. If the operations is not in accordance with the instructions or without the company's authorization units and individuals to open are not included in this scope.

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